


ScanAfia X ICU hospital bed

Versatile ICU-bed with extensive features

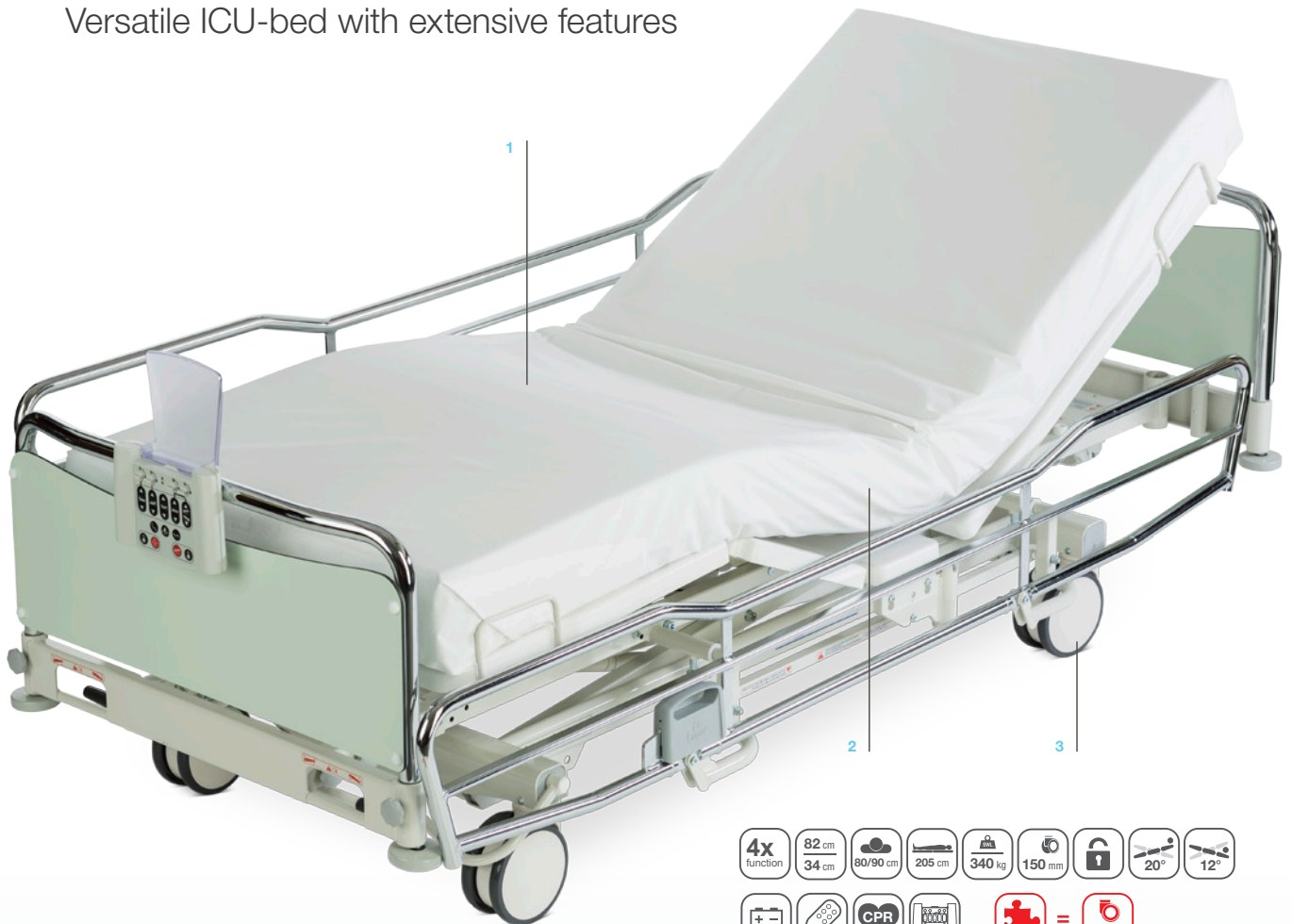

Main features:

1. Very high load capacity: SWL 340 kg
2. Very easy to manoeuvre thanks to 150 mm twin castors
3. Up to 20° anti-Trendelenburg adjustment for effective heart position

High quality ScanAfia X ICU hospital bed has been specifically designed for use in intensive care units. Practical design features provide excellent price-quality ratio. ScanAfia X ICU hospital bed has been tested to fulfil the requirements of the EN 60601-2-52 standard for medical beds.

Despite its impressive safe working load (SWL) of 340kg, the bed is light and its standard 150mm twin castors make it easy to move. The bed is designed to support ergonomic working and patient comfort and accessibility: its exceptionally low minimum height makes work easier for nurses, but also makes getting in and out of bed easier for patients. The back section slide mechanism prevents the patient from sliding towards the foot end of the bed and thus improving patient comfort. A nurse control panel allows bed functions to be locked off from patients. The bed includes special electrical functions such as Trendelenburg, CPR, and cardiac positions.

LOJER® For easy care


Technical data and features	ScanAfia X ICU
Total width (with siderails)	89,5 cm / 99,5 cm
Patient surface width	80 cm / 90 cm
Total length	215 cm
Patient surface length	205 cm
Height adjustment	34 - 82 cm
Recommended mattress size	205 x 78/88 x 13 cm
Total weight (without rails)	89,5 cm: 109 kg 99,5 cm: 111 kg
Castors	Ø150 mm central locking twin castors
Directional wheel	Yes
Safe Working Load (SWL)	340 kg
Electrical system	24V
Back section adjustment	70°
Thigh section adjustment	30°
Calf section adjustment	11°
Trendelenburg / anti-Trendelengburg	12° / 20° (electric)
Removable bed ends	White, blue, yellow and light green
Side rails	Chrome plated or painted: normal and long versions
Other features	Battery backup, hand control, nurse control panel, back section slide, electrical CPR-function, manual emergency lowering of back section, fixing hooks for supporting belts, roller bumpers
CE marked	Yes
Warranty	2 years
Country of Manufacture	Finland (Lojer Oy)

Factory installed options:

- 5th wheel to assist moving

Other accessories:

- Special length (extendable +15 cm) model
- Extra height side rails
- Lifting pole with triangle handle
- Lift support
- IV-pole with one hand quick adjustment
- Auxiliary tray (for linen)
- Bed end (head section)
- Bed end with monitor tray
- Urine bag holder
- Reading and dining tray, on top of the side rails
- Flexible hand controller holder

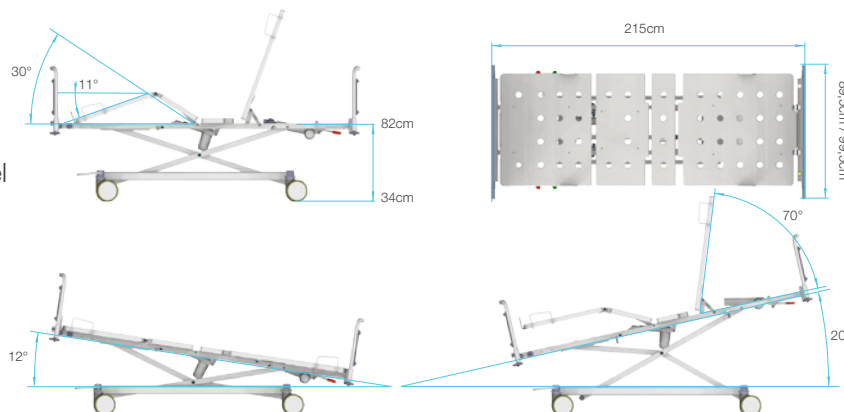

LOJER® For easy care